

U.S. DEPARTMENT OF COMMERCE  
National Technical Information Service

AD-A033 329

## NITROGEN TRIFLUORIDE ANALYTICAL PROCEDURES

AIR FORCE ROCKET PROPULSION LABORATORY  
EDWARDS AIR FORCE BASE, CALIFORNIA

DECEMBER 1976

Reproduced From  
Best Available Copy

ADA033329

ADP-72-18-00

# INTERNAL SECURITY PROTECTIVE PROCEDURES

## FINAL REPORT

AUTHORS: LOUIS A. BEE  
WILLIAM T. LEVOCH, LT. USAF  
DECEMBER 1976

APPROVED FOR PUBLIC RELEASE  
DISTRIBUTION UNLIMITED

AIR FORCE ROCKET PROPULSION LABORATORY


DIRECTOR OF SCIENCE AND TECHNOLOGY

AIR FORCE SYSTEMS COMMAND

WRIGHT AFB, CALIFORNIA 92513

REPRODUCED BY  
NATIONAL TECHNICAL  
INFORMATION SERVICE  
U.S. DEPARTMENT OF COMMERCE  
SPRINGFIELD, VA 22161

Reproduced From  
Best Available Copy

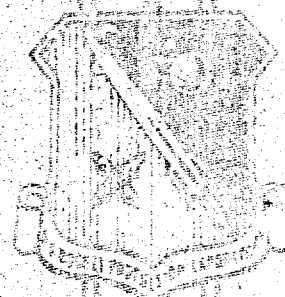

## NOTICES

When U.S. Government drawings, specifications, or other data are used for any purpose other than a definitely related government procurement operation, the Government thereby incurs no responsibility nor any obligation whatsoever, and the fact that the Government may have formulated, furnished, or in any way supplied the said drawings, specifications, or other data, is not to be regarded by implication or otherwise, in any manner licensing the holder or any other person, corporation, or company any rights or permission to manufacture, use, or sell any patented invention that may be related thereto.

## FOREWORD


This report was prepared by the Air Force Rocket Propulsion Laboratory, Edwards AFB, California under Job Order No. 335810S0.

This report has been reviewed by the Information Office/DOZ and is releaseable to the National Technical Information Service (NTIS). At NTIS it will be available to the general public, including foreign nations.

This report is unclassified and suitable for general public release.



LOUIS A. DEE  
Author

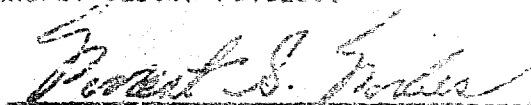

FORREST S. FORBES, Chief  
Propellants Section



WILLIAM T. LEYDEN, 1st Lt, USAF  
Author



EDWARD E. STEIN  
Deputy Chief, Liquid Rocket Division

Reproduced From  
Best Available Copy

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

| REPORT DOCUMENTATION PAGE  |                       | READ INSTRUCTIONS<br>BEFORE COMPLETING FORM  |
|--|-----------------------|--|
| 1. REPORT NUMBER<br>AFRPL-TR-76-89   | 2. GOVT ACCESSION NO. | 3. RECIPIENT'S CATALOG NUMBER  |
| 4. TITLE (and Subtitle)<br><br>NITROGEN TRIFLUORIDE ANALYTICAL PROCEDURES  |                       | 5. TYPE OF REPORT & PERIOD COVERED<br>Final Report                                 |
|  |                       | 6. PERFORMING ORG. REPORT NUMBER   |
| 7. AUTHOR(s)<br>Louis A. Dee<br>William T. Leyden, III, Lt, USAF   |                       | 8. CONTRACT OR GRANT NUMBER(s)   |
| 9. PERFORMING ORGANIZATION NAME AND ADDRESS<br>Air Force Rocket Propulsion Laboratory<br>Edwards AFB, California 93523   |                       | 10. PROGRAM ELEMENT, PROJECT, TASK<br>AREA & WORK UNIT NUMBERS<br><br>JON 305810SQ |
| 11. CONTROLLING OFFICE NAME AND ADDRESS  |                       | 12. REPORT DATE<br>December 1976   |
|  |                       | 13. NUMBER OF PAGES<br>23  |
| 14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)  |                       | 15. SECURITY CLASS. (of this report)<br><br>UNCLASSIFIED                           |
|  |                       | 15a. DECLASSIFICATION/DOWNGRADING<br>SCHEDULE                                      |
| 16. DISTRIBUTION STATEMENT (of this Report)<br>APPROVED FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED   |                       |  |
| 17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)   |                       |  |
| 18. SUPPLEMENTARY NOTES  |                       |  |
| 19. KEY WORDS (Continue on reverse side if necessary and identify by block number)<br>Nitrogen trifluoride<br>analysis<br>Trace impurities<br>Gas Chromatography   |                       |  |
| 20. ABSTRACT (Continue on reverse side if necessary and identify by block number)<br>A gas chromatograph technique for the analysis of $\text{NF}_3$ is described wherein $\text{N}_2$ , $\text{CO}(\text{O}_2)$ , $\text{NO}$ , $\text{CF}_4$ , $\text{NF}_3$ , $\text{CO}_2$ , $\text{N}_2\text{O}$ , $\text{SF}_6$ and possibly $\text{N}_2\text{F}_2$ are separated isothermally on a single column. In addition, a simple total reactive fluoride measurement is presented. |                       |  |

# TABLE OF CONTENTS

| <u>SECTION</u> |                                  | <u>Page</u> |
|----------------|----------------------------------|-------------|
| I              | INTRODUCTION . . . . .           | 1           |
| II             | EXPERIMENTAL PROCEDURE . . . . . | 3           |
| III            | RESULTS AND DISCUSSION . . . . . | 8           |
| IV             | CONCLUSION . . . . .             | 13          |
| V              | LITERATURE CITED . . . . .       | 19          |

## LIST OF ILLUSTRATIONS

| <u>FIGURE</u> |   | <u>Page</u> |
|---------------|---|-------------|
| 1             | Gas Transfer Manifold . . . . .                       | 4           |
| 2             | Total Active Fluoride Apparatus . . . . .             | 5           |
| 3             | Trace Gas Mixture/Helium Matrix . . . . .             | 9           |
| 4             | NF <sub>3</sub> . . . . .                             | 10          |
| 5             | Gas Chromatogram of Calibration Gas Mixture . . . . . | 12          |
| 6             | Gas Chromatogram of NF <sub>3</sub> . . . . .         | 13          |
| 7             | HF Absorbance <u>vs</u> Concentration . . . . .       | 15          |

|                                 |               |                                     |
|---------------------------------|---------------|-------------------------------------|
| ACCESSION BY                    |               |                                     |
| DATE                            | DATE RECEIVED | <input checked="" type="checkbox"/> |
| BY                              | DATE RECEIVED | <input type="checkbox"/>            |
| CLASSIFICATION                  |               |                                     |
| AUTHORITY                       |               |                                     |
| BY                              |               |                                     |
| DISTRIBUTION/AVAILABILITY CODES |               |                                     |
| NO.                             | AVAIL. CODE   | SPECIAL                             |
| A                               |               |                                     |

## INTRODUCTION

Interest in the use of nitrogen trifluoride ( $\text{NF}_3$ ) as a fluorine oxidizer has recently dictated the need for a military propellant specification. Since the specification represents the portion of a procurement contract which defines the quality of the product, it is necessary that the quality tests be reliable and, at the same time, as simple as possible so that a high quality product is assured without unnecessarily high analytical costs.

$\text{NF}_3$  is a relatively inert compound at ambient temperature, however, the possible impurities can vary from highly reactive (e.g.  $\text{HF}$ ,  $\text{F}_2$ ,  $\text{N}_2\text{F}_4$ ,  $\text{COF}_2$ ,  $\text{N}_2\text{F}_2$ ) to very inert (e.g.,  $\text{N}_2$ ,  $\text{O}_2$ ,  $\text{CF}_4$ ,  $\text{CO}_2$ ,  $\text{N}_2\text{O}$ ). Also impurities of intermediate reactivity (e.g.,  $\text{NO}$ ,  $\text{NO}_2$ ,  $\text{CO}$ ) can be present.

Analysis of such a potentially complex mixture requires some compromise with respect to specificity in order to maintain reasonable analysis costs. A brief inspection of the objectives for the analysis of the propellant reveals that the impurities can be grouped into several categories; (a) those that lower performance by dilution, (b) those that lower performance by interference with the fuel/oxidizer reaction, and (c) those that affect the storability of the propellant. A number of these potential impurities are categorized in Table I.

Table I, Impurities Categorized by Detrimental Effect

| a (diluent)       | b (interference) | c (corrosive)                 |
|-------------------|------------------|-------------------------------|
| N <sub>2</sub>    | O <sub>2</sub>   | HF                            |
| CF <sub>4</sub>   | CO               | F <sub>2</sub>                |
| SF <sub>6</sub> ? | NO               | N <sub>2</sub> F <sub>2</sub> |
|                   | CO <sub>2</sub>  | N <sub>2</sub> F <sub>4</sub> |
|                   | N <sub>2</sub> O | COF <sub>2</sub>              |
|                   | HF               |                               |
|                   | COF <sub>2</sub> |                               |
|                   | NO <sub>2</sub>  |                               |

Group (c) can be determined without specificity since the common active atom is fluorine. A total active fluoride determination will provide a greater penalty for N<sub>2</sub>F<sub>4</sub> contamination than for HF on a mole basis. This imbalance is justified since corrosion potential should be proportional to active fluoride concentration. Total oxidizing capacity of NF<sub>3</sub> impurities has been measured with aqueous I<sup>-</sup>(1) and solid Cl<sup>-</sup>(2); however, in both cases HF was determined separately.

The following describes a less complicated test method for NF<sub>3</sub> using only gas chromatography and a total active fluoride measurement.

## EXPERIMENTAL

Gas transfer manifold. The manifold shown in Figure 1 was fabricated from Hoke type 415 "tee" pattern and Hoke type 413 "straight" pattern valves. The differential pressure gauges are 0 - 760 mm and 0 - 20 mm Wallace and Tiernan types FA 145 and FA 111, respectively. Each component was cleaned with trichlorotrifluoroethylene prior to assembly and the assembled manifold was passivated with HF overnight. All gas mixtures were prepared using this apparatus.

I. R. cell. The I. R. cell used was a 100 mm pathlength micro gas cell equipped with BaF<sub>2</sub> windows (Barnes Engineering, 906-0029) and packless valves (Hoke, 4552Q4M). The cell was also cleaned, followed by passivation with HF.

Total active fluoride. The sample reservoir was assembled from a 280 cu. in. aviator's breathing oxygen tank, Hoke type 413 valves, and a 0 - 760 mm absolute pressure gauge (Matheson, 63 - 5501). The bubbler was fabricated from a Fisher-Porter 100 ml aerosol compatibility tube, a "drilled" Swagelok tee, and 1/4 in. O.D. polyethylene tubing which was heat-sealed at one end and re-opened by drilling several 0.005 in. holes into the heat-sealed area. For increased safety, polytetrafluoroethylene tubing, plugged and drilled, would probably have been a better choice. The bubbler containing 40 ml of 0.1N-NaOH and the sample reservoir are pictured in Figure 2. The apparatus is shown in the configuration used for sampling NF<sub>3</sub>. When HF mixtures were



Figure 1, GAS TRANSFER MANIFOLD


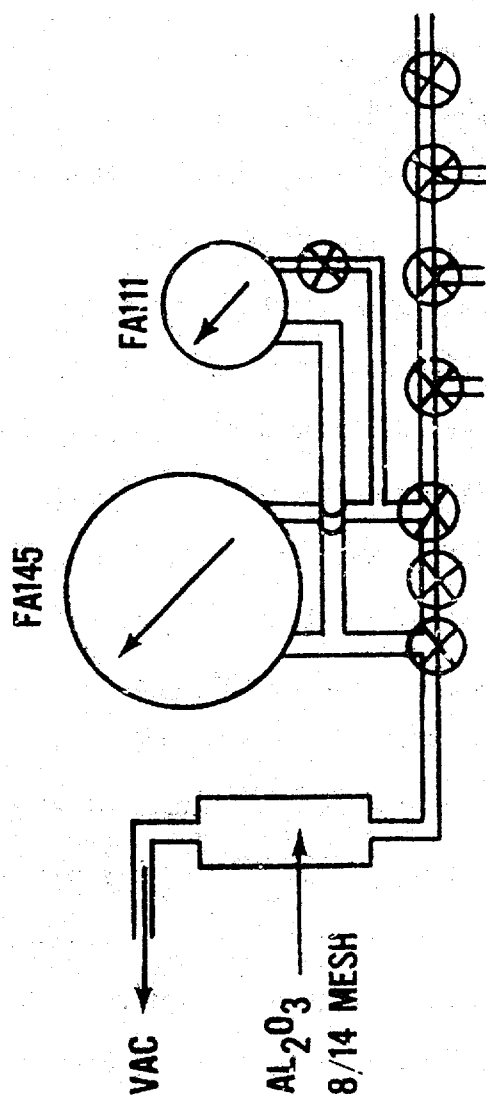


Figure 2, TOTAL ACTIVE FLUORIDE APPARATUS

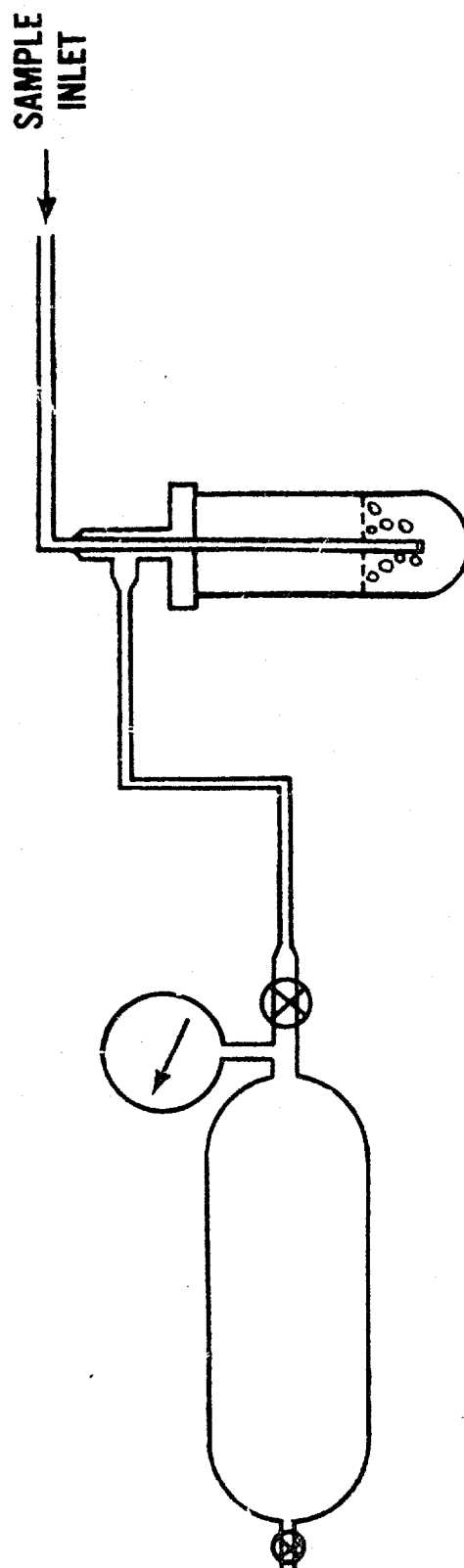

tested, they were prepared in the sample reservoir and drawn through the bubbler, appropriately connected, using a water aspirator. In both cases sample size was computed from the reservoir volume and the pressure change indicated by the gauge. After absorption of the sample, the 0.1N-NaOH was transferred to a vessel containing an equal volume of 0.1M TISAB solution and the fluoride ion concentration was measured using an Orion Fluoride Selective Ion electrode, a Ag/AgCl reference electrode, and an Orion Model 801 Digital pH/MV meter.

Infrared Spectrophotometer. Infrared spectra were obtained using the gas cell described and a Beckman Model 4240 Infrared Spectrophotometer equipped with scale expansion.

Gas Chromatograph. The chromatograms were obtained with a Hewlett-Packard Model 5830A gas chromatograph equipped with a gas sampling valve, a thermal conductivity detector, and a 20 ft. by 1/8 in. O.D. stainless steel column packed with 80/100 mesh Chromosorb 102 (Johns-Manville) or 80/100 mesh Poropak Q (Waters Associates). An active fluoride scrubber was inserted between the sample source and the gas valve. The Scrubber consisted of an 8 in. by 1/4 in. O.D. glass tube packed with 4 in. of 40/50 mesh KI followed by 4 in. of 40/50 mesh  $\text{Na}_2\text{S}_2\text{O}_3 \cdot 5\text{H}_2\text{O}$ . The outlet line of the gas valve was located over a propane/air burner to destroy the excess  $\text{NF}_3$  which was passed through the valve.  $\text{NF}_3$  was analyzed both with and without the active fluoride scrubber in place. Test conditions for the gas chromatograph are listed in Table II.

Table II, Gas Chromatograph Test Conditions

|                     |       |                      |
|---------------------|-------|----------------------|
| Detector type       | ----- | Thermal conductivity |
| Helium flow ml/min. | ----- | 15                   |
| Sample Volume ml    | ----- | 1 @ ambient pressure |
| Column temp. °C.    | ----- | 35                   |
| Inlet Temp. °C      | ----- | 85                   |
| Detector temp. °C   | ----- | 125                  |

Improved resolution between CO and O<sub>2</sub> was observed with Poropak Q at 25°C.

Calibration Standards. All gas mixtures were prepared directly by partial pressure measurement; no serial dilution was used. Trace contaminants were 99.0% pure or better and diluents were greater than 99.99% pure helium or nitrogen.

## RESULTS AND DISCUSSION

Infrared analysis. Figures 3 and 4 are infrared spectra of the calibration gas mixture described in Table III and a sample of  $\text{NF}_3$ , respectively. Both were obtained at 760 mm using the same cell. Data from Pierson, et al (4) and Simons (5) aided in band assignment. Inspection of the figures shows that absorption regions for many of the anticipated contaminants are not located in windows of the  $\text{NF}_3$  spectrum. In addition, poor sensitivity, at the test conditions used, will likely allow only  $\text{NO}_2$ ,  $\text{CF}_4$ , and  $\text{N}_2\text{O}$  to be determined without scale expansion. Although not illustrated by the figures, the useful absorption band for  $\text{N}_2\text{F}_4$ ,  $1270\text{ cm}^{-1}$ (3), is likely overlapped by the  $\text{CF}_4$  band.

Gas chromatograph analysis. A calibration standard containing a variety of compounds anticipated as impurities in  $\text{NF}_3$  was prepared and analyzed using the test conditions described. Table III shows the composition of the mixture and the sensitivities of each component.

Figure 3, TRACE GAS MIXTURE/HELIUM MATRIX (760mm)

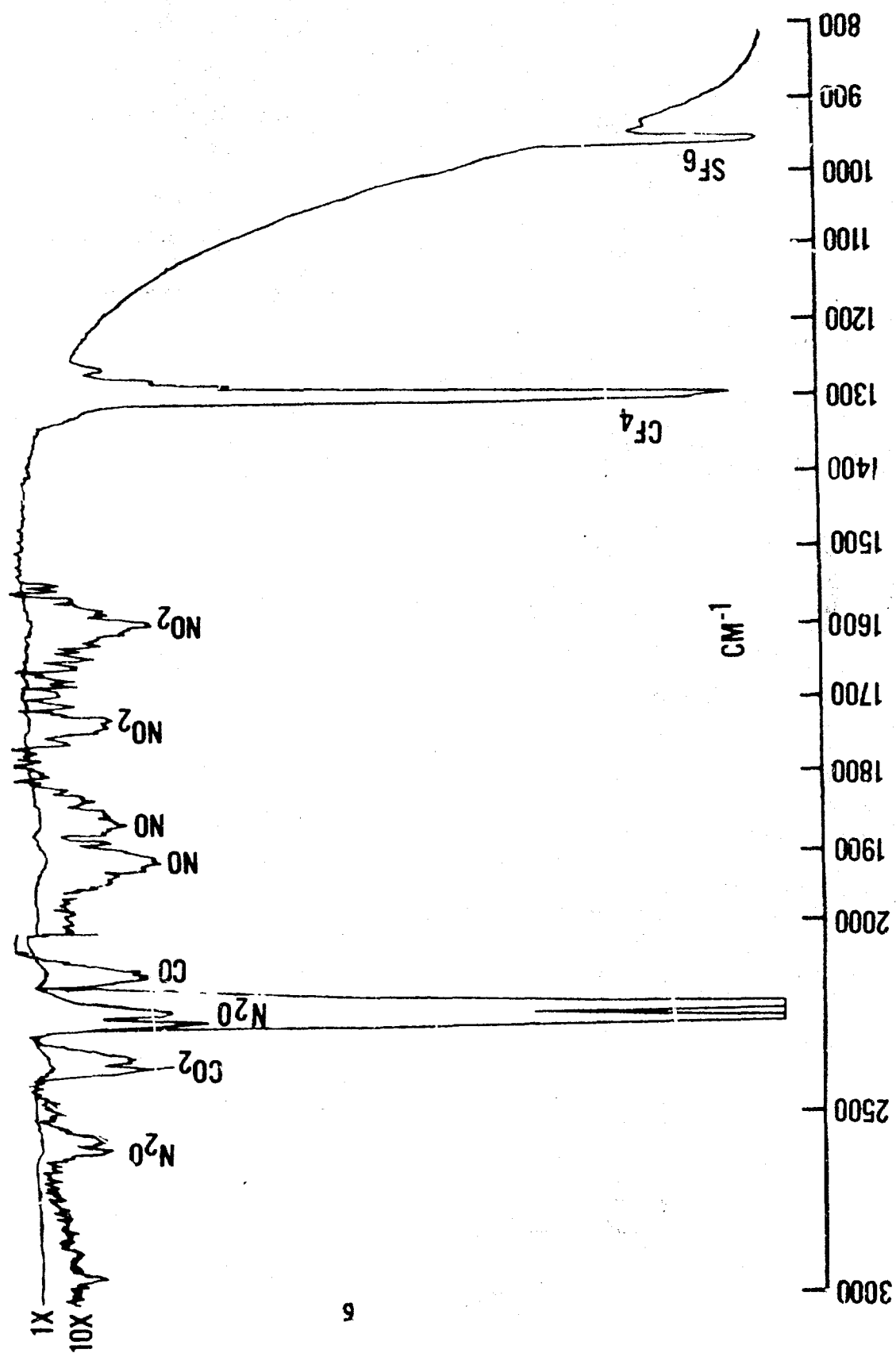


Figure 4, NF<sub>3</sub> (760mm)

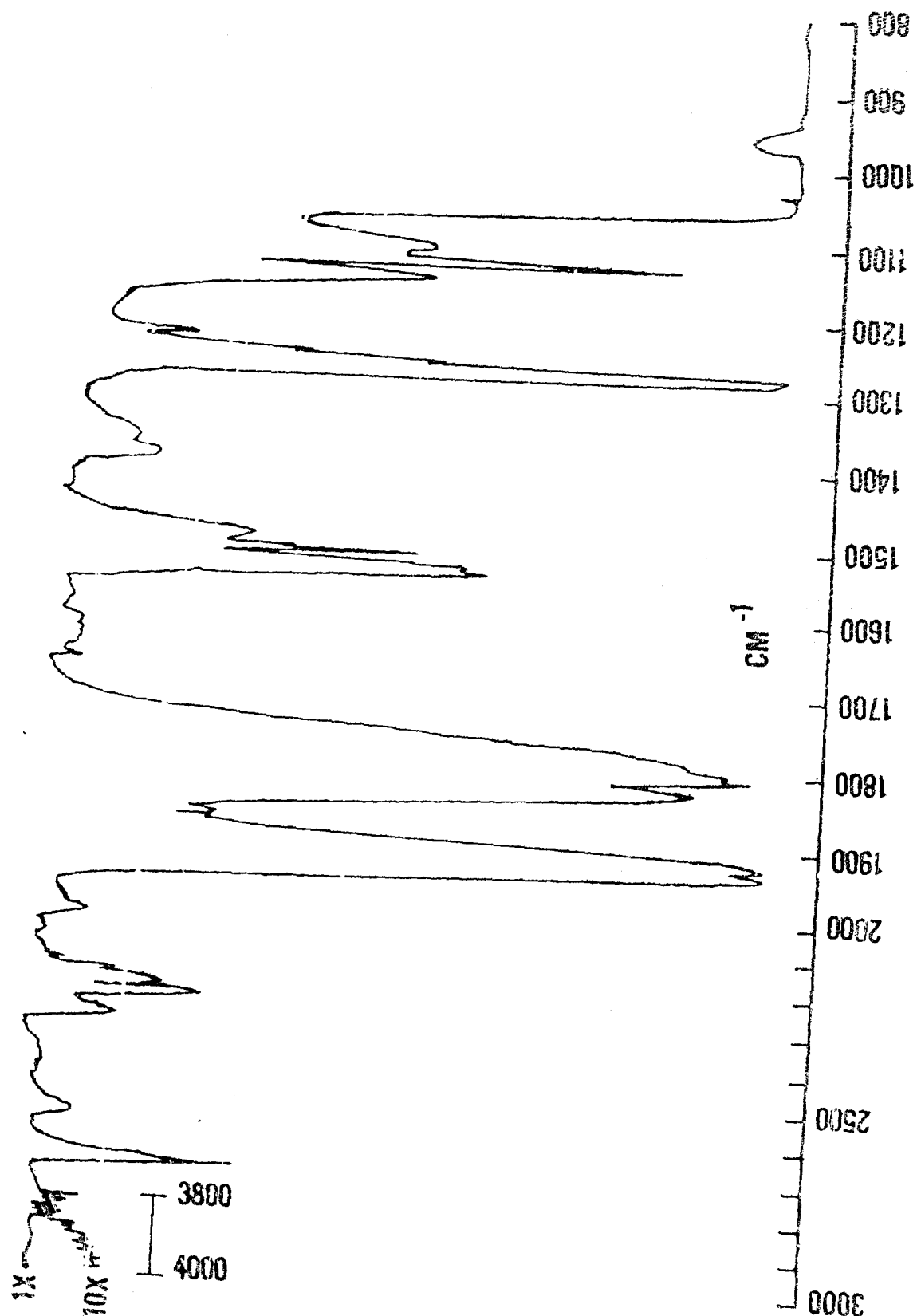


Table III, Calibration Standard Analysis Results

| Component           | Concentration (ppm/vol) | Sensitivity (ppm/area unit) |
|---------------------|-------------------------|-----------------------------|
| N <sub>2</sub>      | 920                     | 0.0710                      |
| CO(O <sub>2</sub> ) | 1010                    | 0.0774                      |
| NO                  | 920                     | 0.0868                      |
| CF <sub>4</sub>     | 920                     | 0.0546                      |
| CO <sub>2</sub>     | 520                     | 0.0597                      |
| N <sub>2</sub> O    | 920                     | 0.0644                      |
| SF <sub>6</sub>     | 920                     | 0.0466                      |

Figure 5 is a typical chromatogram of the calibration mixture, and Figure 6 is a chromatogram of an NF<sub>3</sub> sample using the same test conditions. Cis and trans N<sub>2</sub>F<sub>2</sub> have been reported to appear between the NF<sub>3</sub> and CO<sub>2</sub> peaks (2). Unfortunately, oxygen and carbon monoxide are not well resolved with this column, however, the N<sub>2</sub>, O<sub>2</sub>, and CO<sub>2</sub> peaks can be diverted with a column switching valve to a Molecular Sieve column if oxygen/carbon monoxide separation is desired. The separation is adequate using Poropak Q if the O<sub>2</sub> level is no more than twice the CO level.

Since some of the potential contaminants in NF<sub>3</sub> are strong oxidizers and very corrosive, some means of removing these compounds was desired so that the chromatograph would not be damaged. The KI/Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub> scrubber was installed and the data shown in Table IV was acquired to determine if the scrubber affected the NF<sub>3</sub> analysis results.



FIGURE 5, GAS CHROMATOGRAM OF CALIBRATION GAS MIXTURE

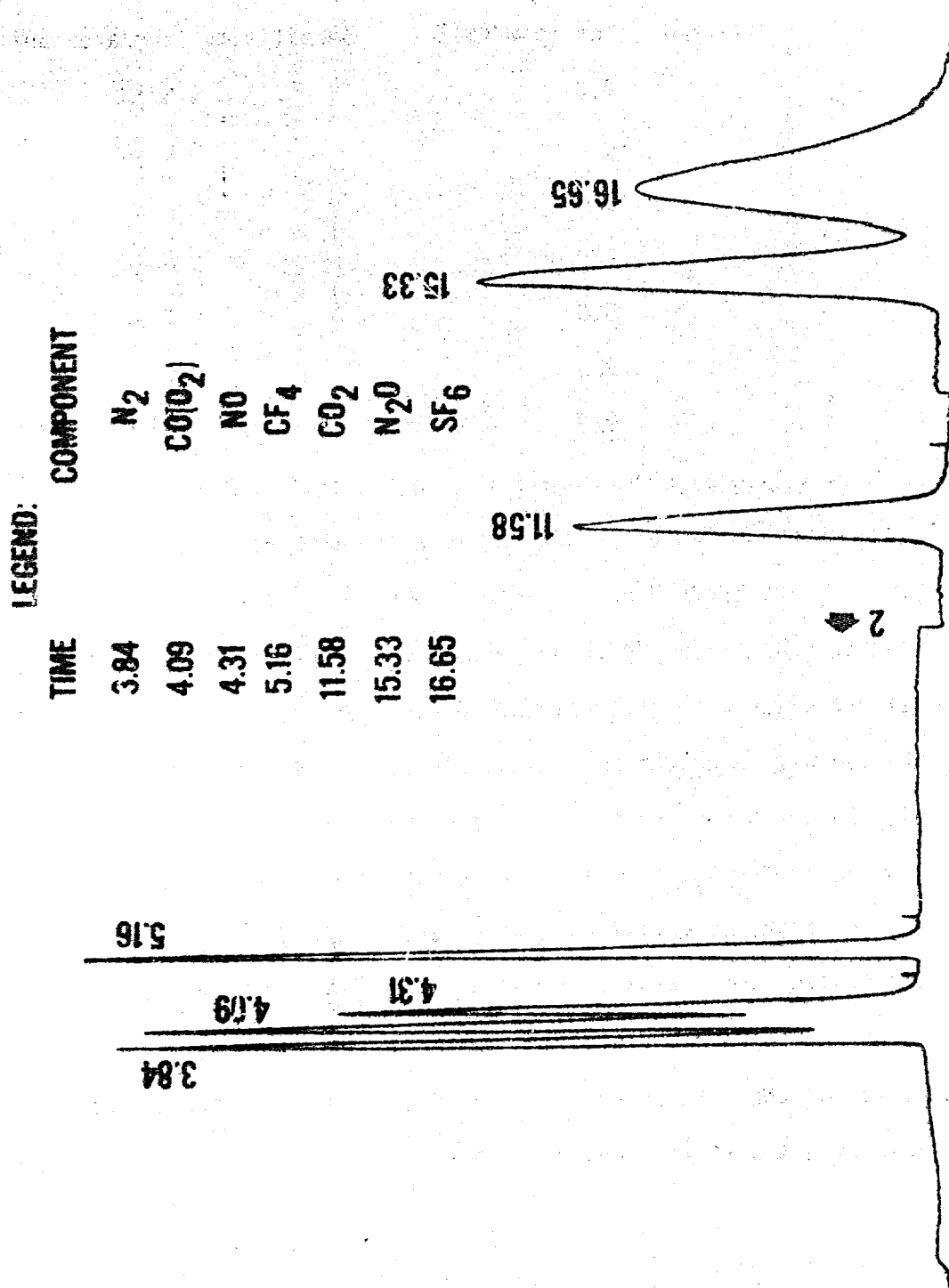


Figure 6, GAS CHROMATOGRAM OF NF<sub>3</sub>

LEGEND:

TIME(MIN.)

COMPONENT

15.18

N<sub>2</sub>

CO(O<sub>2</sub>)

NO

CF<sub>4</sub>

NF<sub>3</sub>

CO<sub>2</sub>

N<sub>2</sub>O

3.80

4.00

4.24

4.99

5.29

11.55

15.18

2

4.99

4.24

4.00

3.80

5.29

11.55

2


Table IV,  $\text{NF}_3$  Analysis Results

| Component                | No Scrubber(%) | Scrubber(%) | Replicate range<br>(6 determinations) |
|--------------------------|----------------|-------------|---------------------------------------|
| $\text{N}_2$             | 0.23           | 0.21        | 0.04                                  |
| $\text{CO} (\text{O}_2)$ | 0.02           | 0.01        | 0.01                                  |
| $\text{NO}$              | 0.02           | 0.02        | 0.00                                  |
| $\text{CF}_4$            | 0.30           | 0.30        | 0.00                                  |
| $\text{NF}_3$            | 99.35          | 99.39       | 0.05                                  |
| $\text{CO}_2$            | 0.007          | 0.008       | 0.001                                 |
| $\text{N}_2\text{O}$     | 0.08           | 0.08        | 0.00                                  |

Some minor differences are indicated, however, they are not significantly different from the replicate variation. If higher levels of oxidizing impurities were present in the  $\text{NF}_3$  sample, use of the scrubber might affect the indicated levels of some measured impurities by release of similar reaction products (i.e.,  $\text{N}_2$  from  $\text{N}_2\text{F}_4$  and  $\text{N}_2\text{F}_2$ ,  $\text{CO}$  from  $\text{COF}_2$ , or  $\text{NO}$  from  $\text{NOF}$ ). Thus,  $\text{NF}_3$  which contains these reactive fluorine compounds may be twice penalized for their presence.

Reactive fluoride analysis. Use of the infrared spectrophotometer for determination of  $\text{HF}$  in  $\text{NF}_3$  is relatively common (1, 2, 6). Figure 7 is a typical calibration curve using data from standards prepared individually by partial pressure measurement in the 100 mm pathlength cell. It is apparent, from the low absorbance values, that scale expansion is useful at the noted test conditions. Other reactive fluorine compounds (i.e.,  $\text{F}_2$ ,  $\text{N}_2\text{F}_4$ ,  $\text{NOF}$ ,  $\text{N}_2\text{F}_2$ ,  $\text{COF}_2$ , etc.) have been determined by combinations of infrared spectrophotometry, red./ox, titrimetry (1) and oxidation of  $\text{Cl}^-$  to  $\text{Cl}_2$  with

Figure 7. HF ABSORBANCE VS. CONCENTRATION

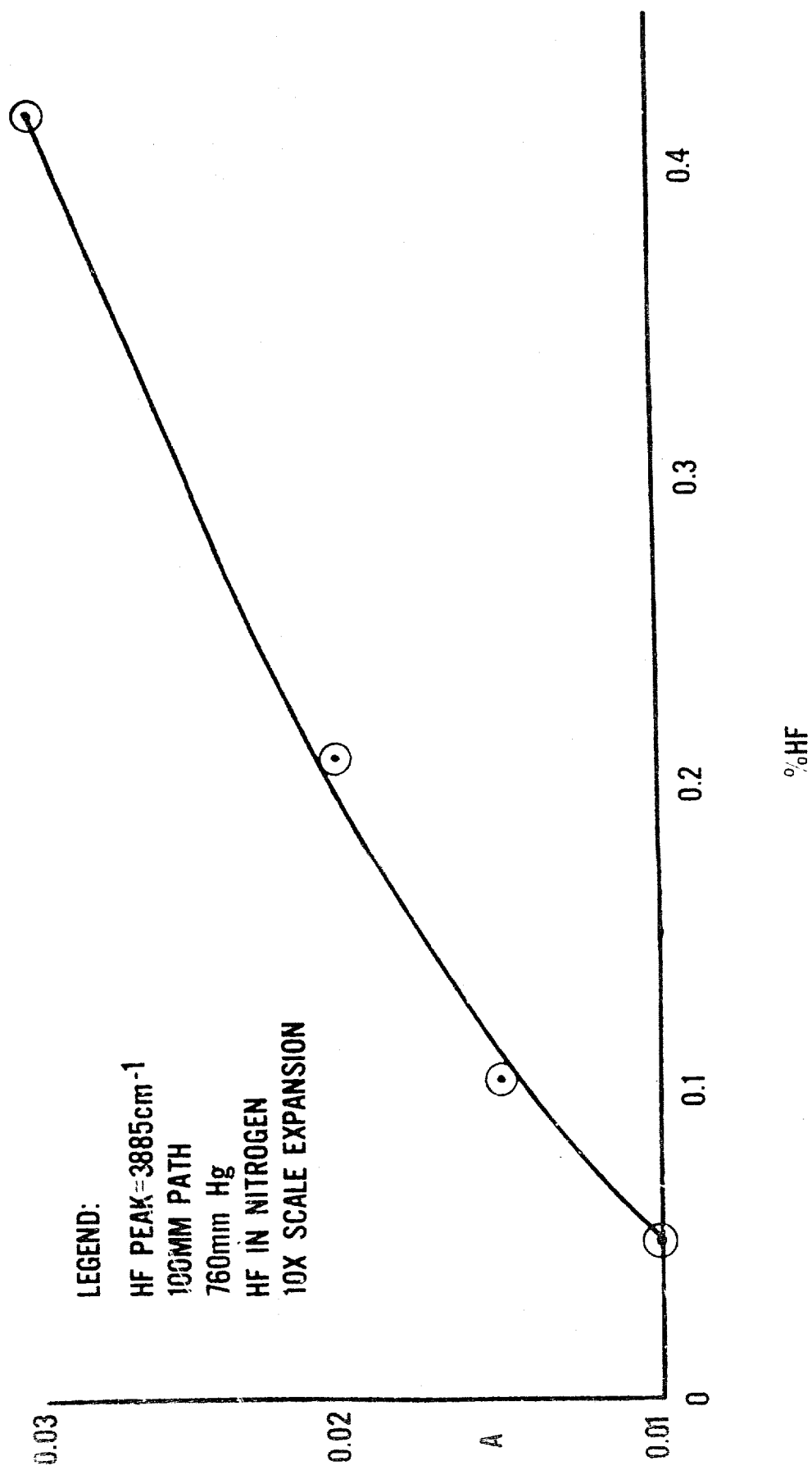

subsequent gas chromatography analysis (2). These techniques are specific for total oxidizing species, however, their use necessitates several separate analyses and increased costs. If the primary reason for controlling these reactive species in  $\text{NF}_3$  is to minimize corrosion then it seems a total reactive fluoride measurement would suffice.

Aqueous 0.1N-NaOH should hydrolyze the reactive fluorine compounds, leaving fluoride ion in proportion to the concentration and molecular formula of each specie. This will result in a graduated penalty for corrosion potential relative to the species present. This appears to be a desirable feature of the measurement.

The apparatus shown in Figure 2 was evaluated for collection efficiency with HF standards in nitrogen,  $\text{NF}_3$ , and HF in  $\text{NF}_3$ . Table V lists the data.

Table V, Evaluation of Total Active Fluoride Apparatus

| Sample            | %HF (calculated) | %HF (found) |
|-------------------|------------------|-------------|
| HF/ $\text{N}_2$  | 0.14             | 0.13        |
| HF/ $\text{N}_2$  | 0.19             | 0.19        |
| HF/ $\text{N}_2$  | 0.14             | 0.14        |
| HF/ $\text{N}_2$  | 0.12             | 0.15        |
| HF/ $\text{N}_2$  | 0.25             | 0.25        |
| $\text{NF}_3$     | -                | 0.004       |
| $\text{NF}_3^*$   | -                | 0.000       |
| HF/ $\text{NF}_3$ | 0.21             | 0.22        |

\*The  $\text{NF}_3$  from the previous test was passed through a second bubbler

Using the test conditions described earlier, the final fluoride ion concentration is ~50ppm/wt. if the gas sample contains 0.1%HF. Table V data was obtained using sample flows through the bubbler of ~100 ml/min., thus indicating that the reaction is rapid and no losses are apparent.

## CONCLUSION

A simple analysis technique for  $\text{NF}_3$  has been demonstrated wherein a relatively complete quantitative description can be obtained with inexpensive laboratory equipment used with reasonable care.

#### LITERATURE CITED

1. Product Specification for  $\text{NF}_3$ , Code Ident. 11982, Draft, TRW Systems Group, Redondo Beach, CA.
2. Acceptance Test Procedure for  $\text{NF}_3$ , No. 7836, North American Rockwell, Rocketdyne Division, Canoga Park, CA.
3. Correspondence, E. A. Burns, Letter No. 4341.75-277, 25 Nov 1975, Applied Chemistry Dept., TRW Systems Group, Redondo Beach, CA.
4. "Fluorine Chemistry", Vol. II, J. H. Simons, ed., 1954, Academic Press Inc., New York, N. Y.
5. R. H. Pierson, A. N. Fletcher, E. S. Gantz, Anal. Chem, 28, 1218 (1956).
6. Mil-P-27405, Propellant, Fluorine, 28 June 1968.